

PRESHAPED TRACHEOSTOMY TUBES

TRACHEOFIX

TRACHEOSTOMY TUBE MADE OF THERMOSENSITIVE PVC

The TracheoFix tube can be used universally in intensive therapy, ENT and long-term care. As a result of the innovative adaptation and connecting mechanism for inner

cannula systems, neither the fixation flange nor the tracheostoma are touched when the tube is exchanged.

- with cuffed versions:
 - valve for Luer and Luer-lock syringe tips
 - blue pilot balloon
 - low-pressure cuff
- with or without phonation window
- multifunctional connector set with rotary locking mechanism (depending on version with phonation valve, cough cap, sealing cap, 15 mm connector, etc.)
- soft fixation flange, flexible in all directions, with letter coding
- fixation flange and inner cannula colour-coded according to size
- innovative connection mechanism for inner cannula systems (rotary locking)
- continuous X-ray marker along the tube shaft
- latex-free
- sterile
- single patient use

ACCESSORIES

Accessories for TracheoFix Set available separately.

See page 15


SPECIAL FEATURE

TracheoFix Set for PDT.

See page 21

SIMPLE AND SAFE IDENTIFICATION


1) Marking letters on the fixation flange:

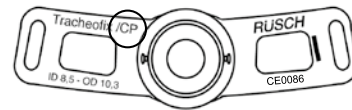


- C = Tube with cuff (without phonation)
- P = Tube with phonation (without cuff)
- CP = Tube with cuff and phonation
- = Tube without cuff and phonation

2) Coloured circle on the inner cannula and colour-code according to size:

interrupted coloured circle:
with phonation

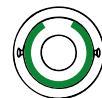

closed coloured circle:
without phonation

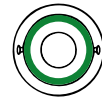

TRACHEOFIX TUBE

RUSCH

ORDER SIZE	I.D.	O.D.	TOTAL LENGTH (CLL=A+B+C)	DIM. A	DIM. B	DIM. C	ANGLE θ (APPROX.)	CUFF Ø
● 7.0	7.2 mm	9.0 mm	70.0 mm	–	–	70.0 mm	102°	27.0 mm
● 8.5	8.7 mm	10.3 mm	78.0 mm	–	–	78.0 mm	94°	27.0 mm
● 10.0	10.2 mm	12.3 mm	86.0 mm	–	–	86.0 mm	90°	32.0 mm
● 11.0	11.2 mm	13.3 mm	86.0 mm	–	–	86.0 mm	92°	35.0 mm

TRACHEOFIX TFC SET
with cuff, without phonation window



- SET COMPONENTS**
- tracheostomy tube
 - 2 colour-coded inner cannulas with 15 mm connector, premounted fixation ring, (long)
 - comfort neck band
 - insertion aid (obturator)

REF.	ORDER SIZE/I.D.	QTY
121400	7.0 – 11.0	5

TRACHEOFIX TFCP SET
with cuff and phonation window



- SET COMPONENTS**
- tracheostomy tube
 - 4 colour-coded inner cannulas
 - 1 inner cannula with 15 mm connector, premounted fixation ring, (long)
 - 1 inner cannula with low-profile connector, (short)
 - 1 inner cannula with 15 mm connector, phonation, premounted fixation ring, (long)
 - 1 inner cannula with low-profile connector, phonation, (short)
 - fixation ring
 - cough cap
 - phonation valve with O₂-adapter
 - comfort neck band
 - insertion aid (obturator)

REF.	ORDER SIZE/I.D.	QTY
121410	7.0 – 11.0	5

TRACHEOFIX TFP SET
without cuff, with phonation window



- SET COMPONENTS**
- tracheostomy tube
 - 4 colour-coded inner cannulas
 - 1 inner cannula with 15 mm connector, premounted fixation ring, (long)
 - 1 inner cannula with low-profile connector, (short)
 - 1 inner cannula with 15 mm connector, phonation, premounted fixation ring, (long)
 - 1 inner cannula with low-profile connector, phonation, (short)
 - fixation ring
 - cough cap
 - phonation valve with O₂-adapter
 - comfort neck band
 - insertion aid (obturator)

REF.	ORDER SIZE/I.D.	QTY
121420	7.0 – 11.0	5

TRACHEOFIX TF SET
without cuff and phonation window







- SET COMPONENTS**
- tracheostomy tube
 - 2 colour-coded inner cannulas
 - 1 inner cannula with 15 mm connector, premounted fixation ring, (long)
 - 1 inner cannula with low-profile connector, (short)
 - fixation ring
 - cough cap
 - comfort neck band
 - insertion aid (obturator)

REF.	ORDER SIZE/I.D.	QTY
121415	7.0 – 11.0	5





TRACHEOFIX INNER CANNULA

- latex-free
- sterile
- single use

TRACHEOFIX INNER CANNULA						<i>RÜSCH</i>	
	REF.	ORDER SIZE	DESCRIPTION	I.D.	O.D.	QTY	
	121425	● 7.0	<ul style="list-style-type: none"> • 15 mm connector • premounted fixation ring 	5.5 mm	6.7 mm	10	
		● 8.5		7.0 mm	8.2 mm	10	
		● 10.0		8.5 mm	9.7 mm	10	
		● 11.0		9.5 mm	10.7 mm	10	
	121430	● 7.0	<ul style="list-style-type: none"> • low-profile connector • fixation ring 	5.5 mm	6.7 mm	10	
		● 8.5		7.0 mm	8.2 mm	10	
		● 10.0		8.5 mm	9.7 mm	10	
		● 11.0		9.5 mm	10.7 mm	10	
	121435	● 7.0	<ul style="list-style-type: none"> • 15 mm connector • premounted fixation ring • phonation 	5.5 mm	6.7 mm	10	
		● 8.5		7.0 mm	8.2 mm	10	
		● 10.0		8.5 mm	9.7 mm	10	
		● 11.0		9.5 mm	10.7 mm	10	
	121440	● 7.0	<ul style="list-style-type: none"> • low-profile connector • fixation ring • phonation 	5.5 mm	6.7 mm	10	
		● 8.5		7.0 mm	8.2 mm	10	
		● 10.0		8.5 mm	9.7 mm	10	
		● 11.0		9.5 mm	10.7 mm	10	

TRACHEOFIX ACCESSORIES

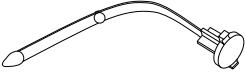- single patient use

TRACHEOFIX ACCESSORIES				
	REF.	ORDER SIZE	DESCRIPTION	QTY
	121475	7.0	insertion aid	1
		8.5		1
		10.0		1
		11.0		1
	121445	–	sealing cap for • weaning	4
	121450	–	phonation valve • for 15 mm connectors with SoftMembrane technology • with O ₂ -adaptor	2
	121455	–	phonation valve • with SoftMembrane technology • with O ₂ -adaptor	2
	121460	–	cough cap	4
	121465	–	fixation ring to secure the TracheoFix inner cannula to low-profile connector	4